

STRENGTH IN MOTION

Designed to MOVE with you.

The QUICKIE® 7 Series redefines ultra lightweight wheelchair technology for the active user. Its customized ergonomics, premium materials, and open back design offer you ultimate freedom of movement. Where will you go?




7
SERIES



TOP OF THE LINE MATERIALS AND TECHNOLOGY

The genius behind the 7 Series' ultra lightweight frame is the 7000 Series aerospace aluminum that has been heat treated using ShapeLoc Technology. This process increases strength and allows for thinner, lighter weight tubing. With an 18% increase in the strength-to-weight ratio over previous designs, a 10% increase in stiffness, and a **transport weight starting at 9.9 pounds***, the 7 Series chairs are light, but mighty.

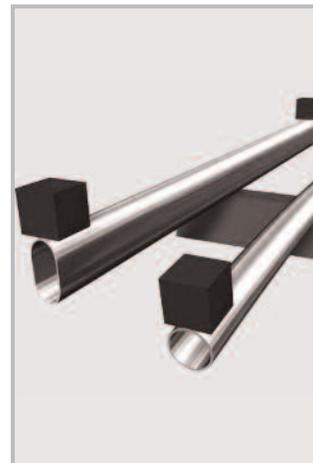

HAND BUILT JUST FOR YOU

As an extension of you, your wheelchair should be custom built for you, your personality, and your lifestyle. The 7 Series offers our highest level of customization with 3 caster settings, 2 frame lengths (per seat depth), 2 frame models, 7 rear wheels, 9 rear tires, and premium selections including the Freestyle Backrest System, Frame Graphics, 4-Link Suspension, QUICKIE Xtender Power Assist, and more.

*Transport weight for the 7RS



Freestyle™ Backrest System

The 7 Series' minimal, ultra-modern design takes backstage so YOU are in the limelight. The Freestyle Backrest System's unmistakable style makes your 7 Series chair a true original. It also breaks barriers to independence by providing a 360° range of upper body motion and lowering the overall weight of your chair.



7^R and 7^{RS} Models

Choose the 7^R if you want to make quick back angle adjustments to tailor your chair to various activities during the day. Choose the 7^{RS} if you know exactly what you need. Its fully rigid backrest and axle plate reduce the overall chair weight by an average of 1 lb. and transfer the energy directly from your arms to the axles for the most efficient propulsion.

CUSTOMIZED
ERGONOMICS,
PREMIUM
MATERIALS,
AND OPEN BACK
DESIGN
FOR THE ULTIMATE
FREEDOM
OF MOVEMENT.



Frame Graphics

Choose from five different QUICKIE exclusive graphics to add style to any frame color, including the 7 Series' exclusive Glossy White. Visit the 7 Series product page on www.SunriseMedical.com to see all graphics available.



Carbon Fiber Accents

Reduce weight and add style with the 7 Series' carbon fiber components. Choose from the carbon fiber footrest, side guards, JAY® J3™ Carbon Fiber Back, or Topolino Carbon Core wheels.



When you want to really break the mold, our Built-4-Me team will help you **Style Your Ride** and create a truly original design with the head-turning style you've imagined.

- Custom wheels
- Colored handrims, accents, and Spinerdy® spokes
- Custom frame colors, upholstery, and JAY® seating



7 SERIES SPECS

HCPCS Code: K0005

Weight

7R: 18.89 lbs.

7RS: 17.67 lbs.

Transport Weight

7R: 11.1 lbs.

(w/o wheel locks: 10.06 lbs.)

7RS: 9.88 lbs.

(w/o wheel locks: 8.84 lbs.)

Weight Capacity

265 lbs.

Seat Width

12" to 20"

Seat Depth

12" to 20"

Center of Gravity Adjustment

0" to 4"

Back Height Adjustment

8" to 20"

Front Seat Height

16" to 21"

Rear Seat Height

13" to 20"

Frame Angle

75°, 80°, 85°

Warranty

Limited/Lifetime